

For that Perfect X-Ray
www.alerio.in

ALERIO[®]
Smart4200
MOBILE X-RAY



Economical, Compact, Powerful HF MOBILE X-RAY



Manufacturer, Sales, Service :
IATOME ELECTRIC (I) PVT. LTD.



AERB Approved



Intertek



BIS Certified



Enquiry & Sales

M: +91 97875 05551 | 88700 11990 | P: (0)422-4220264
Toll-Free: 1800-833-8090 | E: sales@alerio.in | W: www.alerio.in

SPECIFICATION


Model	ALERIO Smart 4200
Generator Type	High Frequency DC Output
Output kV	40kV-120kV in steps of 1kV
Output mA	100mA maximum
Output Power	4.2 kW (100kV, 42mA)
mAs Range	0.1-250mAs
Exposure Time	0.01s – 5.0s
Focal Spot	0.6mm / 1.8mm IEC
Leakage Radiation	< 1mGy in one hour / (<100mR/hr)
Collimator	Manual Collimator with LED bulb
Dimensions Footprint	850mm x 690mm (Floor Area)
Overall Weight	< 80kG
User Interface	Digital with APR and Exposure Switch
Mains Voltage Range	180VAC – 260VAC, 15A

ALERIO Smart 4200 is a compact yet powerful Mobile X-Ray Machine for general radiographic procedures. It uses advanced High Frequency DC Output x-ray generation technology. The overall size and weight of the equipment is low, making it ideal for applications with space constraints. The DC output enables good x-ray imaging at lower exposure settings. The equipment can be powered from a standard power outlet. It has an in-built voltage stabiliser to counter mains voltage variations. The equipment has a selectable second focus for finer low dose imaging. **ALERIO** Smart is designed and Made in India by an ISO 13485 certified company. It comes with warranty and support. The unit is BIS and AERB regulatory certified.

Typical Usage Environment

- Clinics, Laboratories & Hospitals
- ER & Patient Bedside X-Ray
- Remote Locations & Home Visits
- Secondary to Fixed X-Ray

40kV-120kV, 100mA
250mAs, 4.2kW
Dual Focus

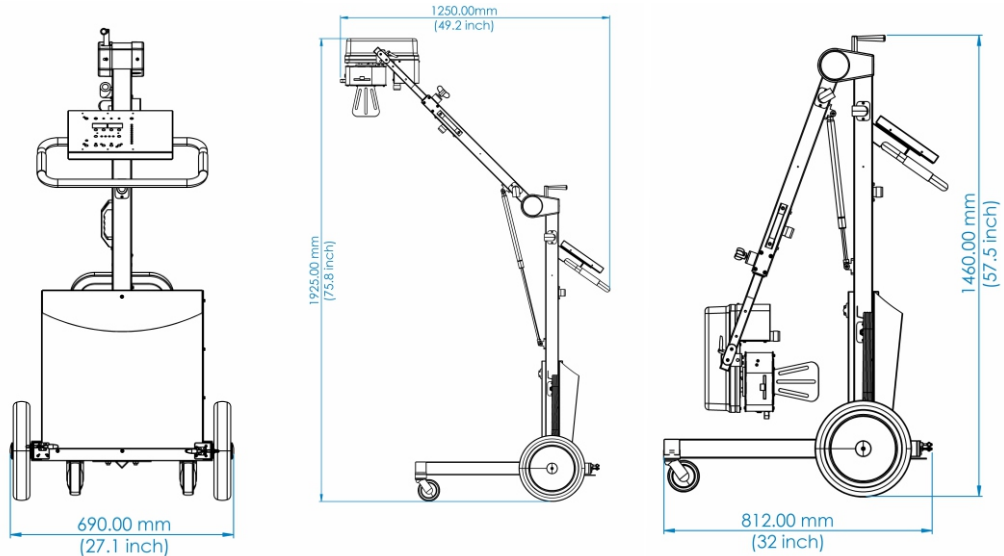

FEATURES

- 100kHz HF DC Generator
- 100mA / 40-120kV / 4.2kW
- Dual Focus LF 1.8mm & SF 0.6mm
- Digital Command & Control
- Compact Size and Weight
- Programmable APR Feature
- In-Built Stabiliser



Images from a ALERIO Smart 4200 mobile x-ray



Manufacturer, Sales, Service :
IATOME ELECTRIC (I) PVT. LTD.
Coimbatore, India - 641 049
www.iatome.in

Enquiry & Sales

M: +91 97875 05551 | 88700 11990 | P: (0)422-4220264
Toll-Free: 1800-833-8090 | E: sales@alerio.in | W: www.alerio.in