

# BR-102 plus

SCHILLER's Ambulatory Blood Pressure Monitoring System combines outstanding accuracy and comfort








## BR-102 plus

# A must-have for ambulatory blood pressure monitoring

The latest version of our established BR-102 plus device results from over 30 years' experience and close collaboration with renowned experts. Ambulatory blood pressure monitoring (ABPM) is considered the gold standard for blood pressure measurement by physicians worldwide. This method provides valuable diagnostic information that stationary hospital or home equipment is unable to provide, such as BP variability, nocturnal BP, as well as nocturnal BP fall and morning BP surge.

SCHILLER's ambulatory blood pressure monitoring system combines two measurement techniques – the auscultatory and the oscillometric – in one: both are used in parallel! This reliable combination ensures correct readings every time.

Furthermore, since the cuff will only inflate as much as necessary, it is hardly possible to hear or feel the BR-102 plus anymore. This added comfort leads to extremely accurate results that can be processed and analysed properly by the evaluation software DARWIN2 ABPM.

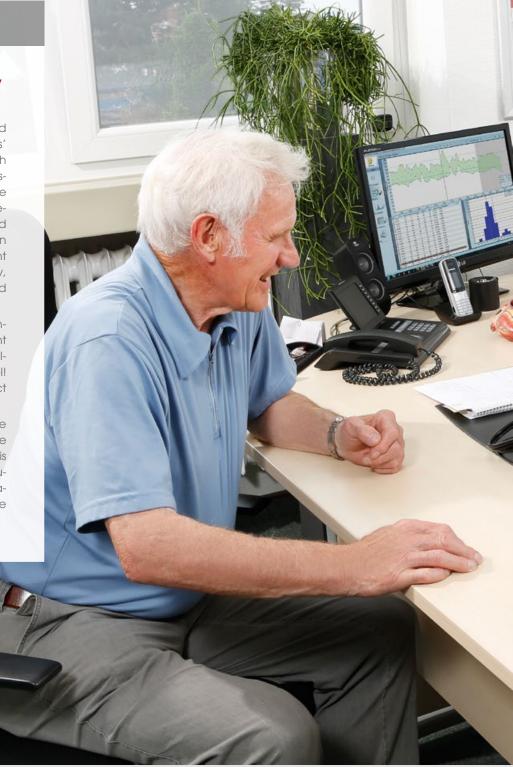

### Accuracy you can rely on

The BR-102 plus is clinically validated and meets all four internationally recognised standards:

- British Hypertension Society (BHS)<sup>1</sup>
- European Society of Hypertension International Protocol (ESH)<sup>1</sup>
- Association for the Advancement of Medical Instrumentation (AAMI) SP10
- Recommended by dabl® Educational Trust
- <sup>1</sup> Test in progress

### Optimal patient comfort for precise results

- Motion-tolerant technology prevents re-inflation and failed readings
- Dynamic inflation to 30 mmHg above the previous systolic reading reduces measurement time
- Built-in pump, and thus almost silent
- Five cuff sizes, from the pediatric to the XL adult size, guarantee optimal measurement accuracy and maximum comfort



## **DARWIN2 ABPM**

# Analysis, interpretation and reporting

DARWIN2 ABPM is a user-friendly, Windows®-based program that allows maximum flexibility for configuration, analysis, interpretation and ABPM reporting. Comprehensive and professional reports can be generated in no time.

### Interpretive summary

Automated analysis of ABPM levels for normality and nocturnal fall according to the AHA $^1$ , ESH $^2$  or JNC7 $^3$  guidelines.

### White-coat analysis

The first-hour data is compared to the full 24-hour results. The software indicates whether the patient presents with a white-coat hypertension.

### One-page report

The one-page report includes ABPM graph, interpretive summary, nocturnal BP fall analysis, and statistics with selected patient data.


### Clear presentation of statistical results

Statistics are provided for awake and asleep phases as well as for all other phases, and include average blood pressure, heart rate, pulse pressure, mean arterial pressure and BP load.

- Recommendations for blood pressure measurement in humans, American Heart Association, 2005.
- $^2\,$  ESH/ESC Guidelines for the management of arterial hypertension, 2013.
- <sup>3</sup> Seventh report of the joint national committee on prevention, detection, evaluation, and treatment of high blood pressure, 2003.

### An ultramodern and versatile device

- Light and compact design
- Simple to set up and operate
- Intended use: white-coat hypertension, drug-resistant hypertension, masked hypertension, hypotensive symptoms with antihypertensive medications

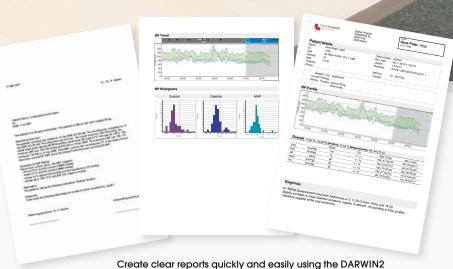

modular reporting system. The level of detail and appearance of reports can be customized at will.





SCHILLER Asia-Pacific / Malaysia

52200 Kuala Lumpur, Malaysia Phone +603 6272 3033 Fax +603 6272 2030 sales@schiller.com.my www.schiller-asia.com



Austria SCHILLER Handelsgesellschaft m.b.H.

A-4040 Linz Phone +43 732 709 90 Fax +43 732 757 000 sales@schiller.at www.schiller.at



### China

Alfred Schiller Medical Equipment Co. Ltd. 100015 Beijing, China Phone +86-10-84565453 info@schillermedical.com www.schillermedical.cn



France SCHILLER Médical S.A.S.

F-67162 Wissembourg/Cedex Phone +33 3 88 63 36 00 Fax +33 3 88 94 12 82 info@schiller.fr www.schiller-medical.com



France (distribution France) SCHILLER France S.A.S. F-77600 Bussy St Georges Phone +33 1 64 66 50 00 Fax +33 1 64 66 50 10 infoschiller@schiller-france.fr www.schiller-france.com



Germany SCHILLER Medizintechnik GmbH

D-85622 Feldkirchen b. München Phone +49 89 62 99 81-0 Fax +49 89 62 99 81-54 info@schillermed.de www.schillermed.de



Hungary SCHILLER Diamed Ltd.

H-1141 Budapest
Phone +36 (1) 383-4780 / 460-9491
Fax +36 (1) 383-4778
sales@schiller.at www.schiller-hungary.hu



### India

SCHILLER Healthcare India Pvt. Ltd. Mumbai - 400 001, India Phone +91 22 6152 3333/ 2920 9141 Fax +91 22 2920 9142 sales@schillerindia.com www.schillerindia.com



### Croatia

Schiller medicinski instrumenti d.o.o.

10000 Zagreb Phone +385 1 309 66 59 Fax +385 1 309 66 60 info@schillerzg.hr www.schiller.ch



### Latin America

SCHILLER Latin America, Inc.

Doral, Florida 33172 Phone +1 954 673 0358 Fax +1 786 845 06 02 info@schillerla.ch www.schillerla.ch



### Poland

SCHILLER Poland Sp. z o.o.

PL-02-729 Warszawa Phone +48 22 8432089 Fax +48 22 8432089 schiller@schiller.pl www.schiller.pl



### Russia & C.I.S.

SCHILLER AG Rep. office 125124 Moscow, Russia 125124 Moscow, Russia Phone +7 (495) 970 11 33 Fax +7 (495) 970 11 33 mail@schiller-ag.com www.schiller-cis.com



### Serbia

SCHILLER d.o.o.

11010 Beograd Phone +381 11 39 79 508 Fax +381 11 39 79 518 info@schiller.rs www.schiller.rs



### Slovenia SCHILLER d.o.o.

2310 Slovenska Bistrica Phone +386 2 843 00 56 Fax +386 2 843 00 57 info@schiller.si www.schiller.si



Spain SCHILLER ESPAÑA, S.A. E-28230-Las Rozas/Madrid Phone +34 91 713 01 76 Fax +34 91 355 79 33 schiller@schiller.es www.schiller.es



### Switzerland

SCHILLER-Reomed AG

CH-8953 Dietikon Phone +41 44 744 30 00 Fax +41 44 740 37 10 sales@schiller-reomed.ch www.schiller-reomed.ch



## Turkey SCHILLER TÜRKIYE

Okmeydani-Sisli – Istanbul Phone +90 212 210 8681 (pbx) Fax +90 212 210 8684 sales@schiller-turkive.com



www.schiller-turkiye.com

SA SCHILLER America Inc.
Doral, Florida 33172
Phone +1 786 845 0620
Fax +1 786 845 06 02 sales@schilleramerica.com www.schilleramerica.com



